



# SLIDE FOLDER

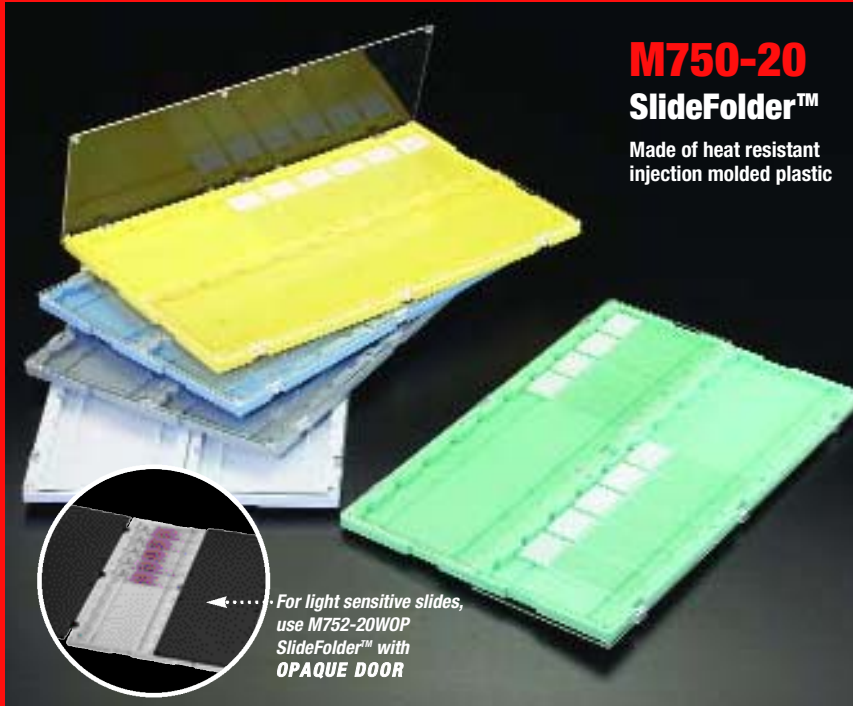

## M750-20 SlideFolder™

Made of heat resistant injection molded plastic

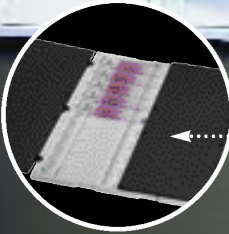

For light sensitive slides, use M752-20WOP SlideFolder™ with OPAQUE DOOR

The SlideFolder™ will hold up to twenty standard microscope slides 75 mm x 25 mm (3 x 1 in.) and is made of two parts: a base holding the slides horizontally offering numbered spaces for easy identification, and transparent doors which can either cover the slides or be swung behind the SlideFolder™ for space saving purposes.


The base is available in 5 different colors, allowing color-coding classification of the slides. Each slide location is identified with a number from 1 to 20. Removal of slides is made easy simply by pressing on one end, which will automatically lift the other end.

The two transparent hinged covers offer a full view of each slide without having to remove it from its position in the SlideFolder™ and allows easy reading of ID labels with or without an optical bar code reader. All units are stackable and take minimum space on laboratory tables or shelves. Will resist temperatures between -80 °C and +100 °C. Not autoclavable.

Dimensions: 192 x 292 x 11 mm H.  
(7 9/16 x 11 1/2 x 7/16 in. H.)

Cat. #	Color	Qty/Cs
M750-20B	Blue	10
M750-20G	Green	10
M750-20GY	Gray	10
M750-20W	White	10
M750-20Y	Yellow	10
M750-20AS	Assorted (two of each color)	10

Cat. #	Color	Qty/Cs
M752-20WOP	Opaque doors and white base	10

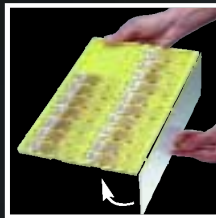

For easier access to slides, simply swing the transparent hinged covers behind the SlideFolder™.



Easy slide removal by pressing down on one end and lifting it from the other.

The SlideTray™ is a convenient microscope slide holder made of heavy-duty plastic lasting many years even under the most adverse conditions. The SlideTray™ will hold up to 20 microscope slides in an almost horizontal position. Each slide can be easily removed and placed back in its position. The SlideTray™ is easily stackable and will take minimum space on any shelf or laboratory counter. Even when trays are stacked, slides are well protected and will not touch the tray above. Available in 4 attractive colors

Dimensions: 206 x 299 x 18 mm H.  
(8 1/8 x 11 3/4 x 11/16 in. H.)

Cat. #	Color	Qty/Cs
M755-20B	Blue	10
M755-20G	Green	10
M755-20W	White	10
M755-20Y	Yellow	10

## M755-20 SlideTray™

Made of heat resistant injection molded plastic



Slides are easily inserted and removed.